



DULET



REF. F360275350

DULET


ULTRALIGHT TACTICAL EVACUATION DEVICE

Ultralight disposable sheet for transfer and coverage, waterproof and absorbent. It keeps the wounded dry to reduce heat loss and protects stretchers and surgical tables to prevent cross-contamination. It can be used to transport, transfer or move casualties of any weight.

FEATURES

- ✓ 100% Latex Free.
- ✓ White upper layer in contact with the patient: 100% hydrophilic polypropylene of 20 g/m².
- ✓ Intermediate layer of S.A.P. sodium polyacrylate average density of 65 g/m² ± 5 g. between two sheets of 5 g cellulose tissue.
- ✓ 42 g/m² green polyester underlay.
- ✓ Approximate weight of the product ready for use 120 g.
- ✓ 50mm 144 decitex tape polyester traction handles sewn to the sheet with 150 dtex black nylon. On this a second one of 25 mm sewn 5 cm leaving sections of 19 cm free.
- ✓ Absorption capacity ≥2000 ml/m².
- ✓ Theoretical traction resistance 220 kg that allows the patient to be moved by pulling the sheet and using it as a transfer.
- ✓ Tolerances: Measurements ± 1 cm/m. Weight ± 10% / pieza.




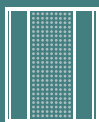

DULET

REF. F360275350

DULET

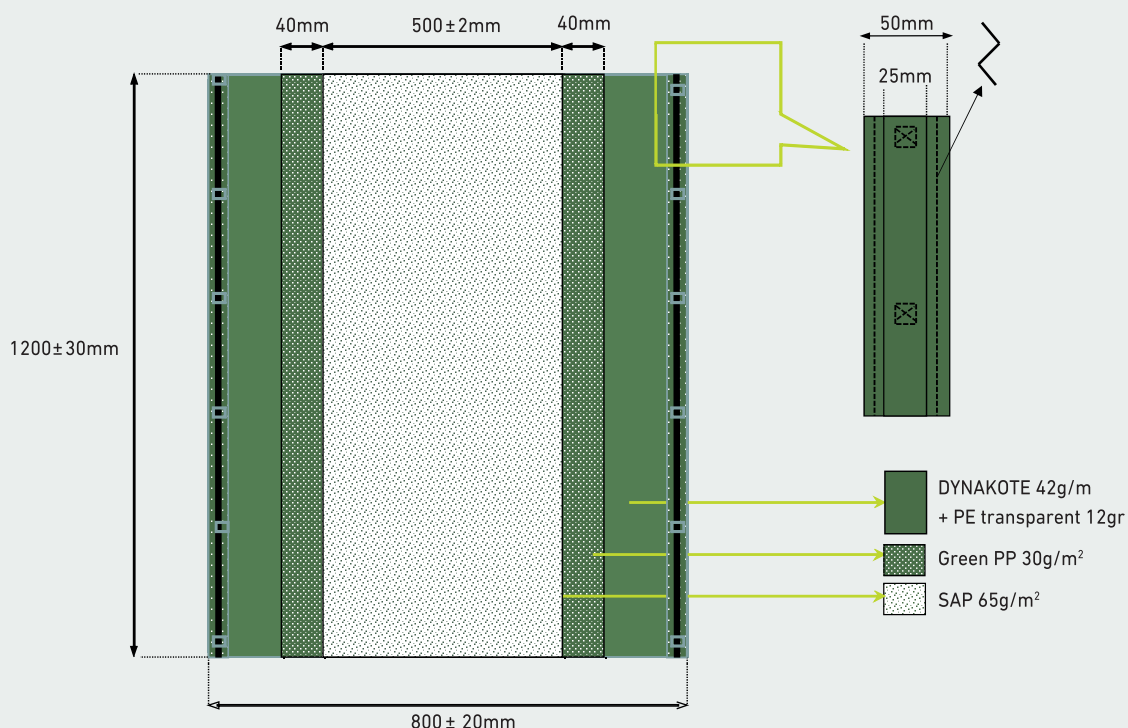

PRESENTATION

120 x 80 cm sheet with 50 x 120 cm absorbent central band.
Individually vacuum packed to carry in a backpack.

CONSERVATION

Once out of its package, keep the product in a dry place.

COMPLIANCE

Non-sterile Class I product. It conforms to EN 13795.
CE self-certification according to European directive 93/42/ECC.

TYPES OF USES

- Rapid extraction of wounded from combat zones.
- Extraction of wounded in cliffs, caves, etc.
- Any situation that requires the rapid transportation of people with disabilities to move on their own.



LATEX FREE



SINGLE USE



WATERPROOF



NON ESTERILE



CE MARKING